




ALCO-SENSOR FST

Evidential Breath Alcohol Tester

Evidential accuracy in the palm of your hand.

The Alco-Sensor FST is approved by the U.S. Dept. of Transportation as an evidential breath testing device and is trusted worldwide by law enforcement agencies as a preliminary breath tester of choice. Easy-to-use and durable, the Alco-Sensor FST provides reliable results for a more accurate, more defensible testing program.

- Accuracy exceeds +/- .005 at .100 BrAC
- Automatic breath sampling feature captures a deep lung breath sample
- Manual breath sampling is available in shy lung situations or when a subject is uncooperative
- Intuitive testing menus guide the operator through a test sequence, reducing operator errors
- Passive testing mode allows you to screen many people quickly for alcohol and is ideal for large events or roadside screening situations
- One-handed operation is valuable for law enforcement officers

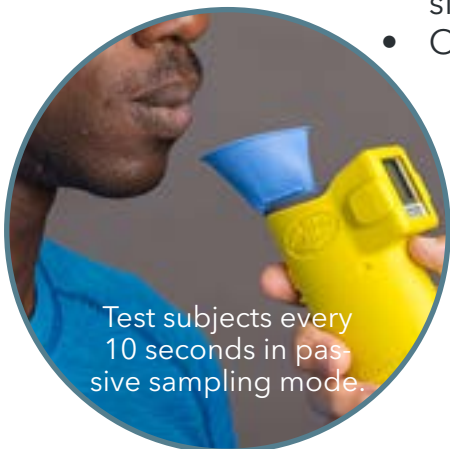

Test subjects every 10 seconds in passive sampling mode.

ALCOPRO
DRUG & ALCOHOL TESTING PRODUCTS

SPECIFICATIONS

Sensor Type:	Electrochemical Oxidation	Calibration:	Microprocessor controlled using gas or wet bath simulator
Measurement Range:	.000 - .040 BrAC	Memory:	Last test
Breath Sampling:	Automatic with manual over-ride	Print Capability:	No
Results Response Time:	5 - 10 sec.	Computer Capability:	No
Power:	2 AA batteries (included)	Operating Temperature:	0 - 50° celcius
Test Modes:	Direct (mouthpiece); Passive (sampling cup)	Size:	5" x 2.5" (6 oz)
Warranty:	One year parts and labor	View our device operation, calibration and accuracy check videos online	

ALCO-SENSOR FST

Item	Standard Product	Intended Markets
207-LE	Alco-Sensor FST, hard-shell case, 25 mouthpieces, battery, sampling cup, drink sniffer, operator manual	Law enforcement, corrections
207-CAN	Alco-Sensor FST, hard-shell case, 25 mouthpieces, battery, sampling cup, drink sniffer, operator manual, mini alcohol gas can	Substance abuse treatment, schools/colleges, non-DOT workplace testing

RELATED ITEMS AND SUPPLIES

Item	Description
343	Alco-Sensor FST mouthpiece
842	TestSafe™ bacterial viral filter
353	Passive sampling cup
324	Alcohol gas tank, 105L
328-06	Regulator 6 LPM



Detects alcohol in beverages.

A Few Good Reasons To Trust Us



Same day shipping
Most orders placed before 2:00 pm EST ship same day.



Knowledgeable Customer Support
Chat online or by phone.



Trusted Supplier more than 40 years.
Supporting our customers since 1982.



Largest inventory
We stock more, so you never run out.