



Premier Biotech — a leading provider of innovative drug testing products and solutions is proud to offer OralTox®, an FDA 510(k) cleared rapid oral fluid testing device. With a continued focus on offering safe and reliable solutions, OralTox delivers accurate and proven results while elevating workplace testing scenarios helping to get candidates to work faster.



### Oral Fluid Testing Advantages

Oral fluid drug testing offers a minimally-invasive, 100% observed collection allowing employers to action immediately with results in just minutes. Using oral fluid, an employer can oversee the candidate as he or she collects their own oral fluid specimens anytime and anywhere thus reducing the likelihood of tampering or cheating the test. The flexibility oral fluid screening provides makes it an ideal solution for onsite drug testing at places of employment, hiring events, and job fairs.


- Expedited hiring with fast test results returned in just minutes
- Less invasive, easy-to-administer procedure that can be done anytime and anywhere
- Observed collection reduces the likelihood of tampering or cheating
- No secured restroom or gender specific requirements

### Why OralTox?

Truly in a league of its own, the streamlined, patented design of OralTox stands out with all the bells and whistles helping to ensure the collection process is easy to administer while highly accurate. OralTox is FDA 510(k) cleared and the only rapid oral fluid device to contain an assay for Oxycodone. When confirmation testing is required (inconclusive results only), simply send the device to the confirmation laboratory. No second collection is required.

- FDA 510(k) cleared device for 8 drugs including, AMP, COC, MTD, MET, OPI, OXY, PCP, THC (Additional panels available)
- Proven accuracy and sensitivity backed by laboratory (LC-MS/MS) results
- Patented design with one drug per strip
- Built-in saturation indicator changes color when a sufficient sample has been collected
- Collector swab meets the International Standard Organization (ISO #10993-1) bio-compatibility standards

### Designed For Speed And Reliability

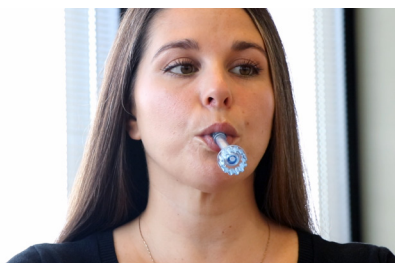

#### Collection

The candidate inserts the OralTox collection swab in their mouth and swabs until the pad is fully saturated. The saturation indicator strip will appear once the proper specimen amount is collected thus removing guess work.



#### Test

With the OralTox device on a flat surface, the collector is inserted with the collection swab face down into the OralTox device. Press down and twist clockwise until the cap of the device is securely locked into place.



#### Results

Within minutes of activating the test, results can be read using the test and control lines on the OralTox strips. If the results are inconclusive, OralTox can be sent to a lab for confirmation testing without the need for a second sample.